

# AVERO Motion E



**beka** DIFFERENT  
BY DESIGN





\* only in combination with a transfer hoist



Promote the mobility and independence of residents with the designer care bathtub AVERO Motion E.

A safe entry is guaranteed by the low entry height.

Partially mobile residents in particular appreciate a relatively independent entry into the care bathtub.



A unique bathing experience for your residents!

Soft shapes, a large amount of space and the ergonomic tub body guarantee a comfortable bathing experience for your residents.

The AVERO Motion E is operated via the practical cable remote control with intuitive symbols.





Sitting position of the tub is ideal for showering your residents



No edges in the door area for unimpeded transfer with the wheelchair



Freely accessible base frame suitable for all modern patient lifts




Large space in the tub body




Allows a clear view of the room without a bulky control panel



Large foot area with a pre-filling volume of 98 liters

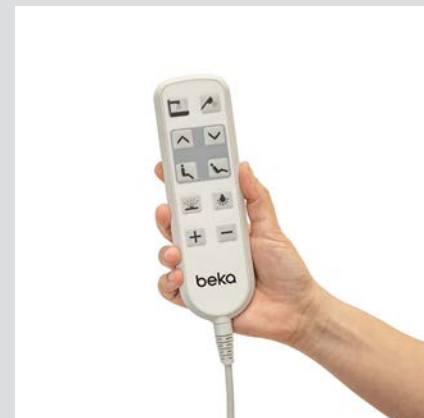

Cable remote control with intuitive symbols for easy operation

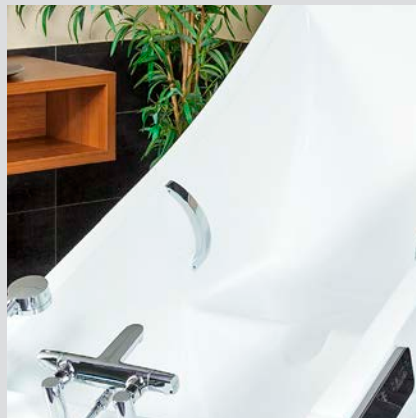

Ergonomic handles provide bathers with optimal support and safety



Special back profile prevents slipping in the tub body



## Elegant design with tilt function

Soft shapes characterize the appearance of the AVERO Motion E. In the sitting position, the tub can be used to shower residents without any problems. The tilt function gently brings residents into the bathing position. Thanks to the large pre-filling volume of 98 l in the foot area of the tub, the residents immediately lie in the warm bath water when the tub tilts back!



## High bathing comfort for every type of resident

The AVERO Motion E offers enough space for tall and heavy residents allowing for a pleasant bathing experience regardless of body size and weight.



## Simple, safe self-entry possible

The large side opening door of the AVERO Motion E enables residents to enter into the bathtub independently. This promotes mobility and gives partially mobile residents a dignified feeling of independence. The low entry height is also ideal for wheelchair transfers.



## Freely accessible base frame

The AVERO Motion E can be freely accessed with all modern bath and patient lifts.






Technical Data


Tub length (entry position / lying position)	1.760 mm / 1.908 mm
Tub width (outside)	860 mm
Lowest entry height	490 mm
Angle of inclination	23°
Tub emptying (liters / min.)	75 l/min
Power connection	230 V, 50 Hz (120 V, 60 Hz)
SWL	350 kg
Pre-filling volume (foot area)	98 l


Product Information


Electrical height adjustment and tub tilting function
Cable remote control with intuitive symbols for all functions
Molded seat with a special back profile to prevent slipping when bathing
Double scalding protection
Ergonomic handles for optimal grip
Seat with hygiene opening
Low entry height of 49 cm


Available Options


- 


Aroma system
- 


Colored-light system
- 

Airspa system
- 

Bluetooth-music system
- 

RAL colors available
- 

Disinfection system
- 

Automatic tub filling
- 

Emergency function in case of power failure



**BEKA Hospitec GmbH** | Am Rübenmorgen 3 | 35582 Wetzlar-Dutenhofen  
Tel +49 641 922 20-0 | E-Mail [info@beka-hospitec.com](mailto:info@beka-hospitec.com)  
[www.beka-hospitec.com](http://www.beka-hospitec.com)

