



Therapy and treatment table MANIPULAT



MANIPULAT 2
Two-part treatment table. Head part positive up to +40° and negative up to -70° adjustable.



MANIPULAT 3 F
Three-part treatment table. Foot part gas-spring supported, positive adjustable by 70°. Middle part is fixed.




MANIPULAT 3 G
Three-part treatment table. Middle and foot part gas-spring supported, positive adjustable by 35°.

Therapy and treatment table **MANIPULAT**

... ideal for demanding treatment methods

The therapy and treatment table MANIPULAT has been designed for treatment methods, which special demands and requirement for high working stress, stability and multi-purpose usage. The height adjustment of the treatment table is absolutely vertical from 480 – 1.000 mm, optional electrical or hydraulic. This offers an ideal adjustment to the therapist and the respective treatment method. The head part can be adjusted positive by +40° and negative by -70° by gas-spring. In case of floor contact of the head part, the gas-spring is automatically released for safety reasons. The treatment area has rounded edges and is padded with fire-resistant artificial leather.

- vertical height adjustment
- high stability and working stress-resistant
- head part positive by +40° and negative by -70° adjustable
- easy patient transfer
- nose-slot integrated in head part
- carrier frame: covered, steel-scissor-system
- all pivotal points teflon-coated

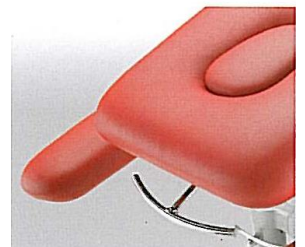

Technical data:

Max. body weight 200 kg

L: 2.050 mm, W: 650 mm, H: 480-1.030 mm

MANIPULAT 2 EL electric height adjustment

Order No. 010721000

MANIPULAT 2 HYD hydraulic height adjustment

Order No. 010720000

MANIPULAT 3 F EL electric height adjustment

Order No. 010723000

MANIPULAT 3 F HYD hydraulic height adjustment

Order No. 010722000

MANIPULAT 3 G EL electric height adjustment

Order No. 010725000

MANIPULAT 3 G HYD hydraulic height adjustment

Order No. 010724000

Power supply for electrical model: 230V, 50Hz

Accessories:

Nose slot

Order No. 010710000

Three-part head part for MANIPULAT

Order No. 010725100

MANIPULAT Castor lifting system

Order No. 010729000

MANIPULAT foot switch

Order No. 010709000

MANIPULAT height adjustment system

Order No. 010734000

Upholstery width 70 cm

Order No. 010727100

Upholstery width 75 cm

Order No. 010727200

Upholstery width 80 cm

Order No. 010727300

Paper holder

Order No. 080701000

Paper roll

Order No. 080702000

Am Rübenmorgen 3
D-35582 Wetzlar-Dutenhofen
Fon +49 (0) 641/92220-0
Fax +49 (0) 641/92220-20
www.beka-hospitec.de

beka
hospitec